

HIWFRS IEC Equipment Inventory

All frontline IEC bags will have a standard stowage. The correct stowage of medical response equipment is detailed in this document.

All medical equipment designated as 'single use' is colour coded **Red**. Equipment that can be cleaned and reused is colour coded **Green**.

To order replacement equipment the following options are available:

1. Equipment marked as **Consumable** can be ordered by emailing '[Station consumable items form](#)' to – stores@hantsfire.gov.uk (**Technical Services**)
2. Equipment marked as '**Indent**' can be replaced by ordering on the Stores Ordering Portal.
3. Equipment marked as '**OET Defect**' should be defected with Control.
4. Replacement **defib pads** are to be ordered via **Technical Services (Stores)**.
5. All **HIWFRS** Medical gas cylinders (**O2/Entonox**) can be exchanged with **Technical Services (Stores)**.
6. **Due to the shelf life of consumable items, stations are advised to use discretion based on operational use/need when maintaining stock levels.**
7. **Consumable items** must be returned to stores **6 months** before expiry date.
8. Blast/OLAES dressings can be stored indefinitely and re-dated as long as the packaging is intact with no loss of vacuum seal.
9. Consumable items marked, **use till expiry date** must be returned to stores (**Aspirin/Glucogel**) or safely disposed of after this date.

Entonox Equipment					
	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Entonox bag	1	11332	OET Defect	N/A
	Entonox cylinder	1	One for one exchange	Control/SCAS (see note 5)	expiry date
	Demand valve & hose assembly	1	11333	OET Defect	N/A
	Filter	5	11359	Consumable	expiry date
	Bite stick	5	11357	Consumable	expiry date

Entonox Equipment cont'd

	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Mask – Child (Size 2)	1	11358	Consumable	Use till expiry date
	Mask – Adult	1	11373	Consumable	Use till expiry date
	Entonox guidance label	1	11644	Consumable	N/A

Infection Control

	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Alcohol hand gel 50ml/350ml	1	11464/10243	Consumable	Use till expiry date
	Bio Guard 500ml Trigger spray	1	11453	Consumable	Use till expiry date
	Clinell Wipes	1	11078	Consumable	Use till expiry date
	FF2 Facemask	1 box 50	12012	Consumable	Dispose after use
	L2 disposable Apron	1	11940	Consumable	Dispose after use

	L3 Disposable gown (with sleeves)	1	12074	Consumable	Dispose after use
	Visor	1	12085	Consumable	Dispose or clean after Use

IEC Bag




	Item	Quantity	Store Number	Replacement	Replacement date
	IEC bag (Total bag comprising three module bags)	1	11336	OET Defect	N/A



	Item	Quantity	Store Number	Replacement	Replacement date
	Module 1 bag	1	11350	OET Defect	N/A
	Module 2 bag	1	11351	OET Defect	N/A
	Module 3 bag	1	11337	OET Defect	N/A
	Double end pouch connector (Co-Responders Only)	1	11612	Consumable	N/A

Module 1 (Diagnostics)

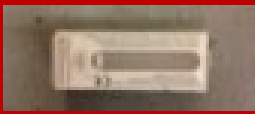

	Item	Quantity	Store Number	Replacement	Replacement date
	BP monitor <i>Netting /Bottom, left side</i>	1	11347	OET Defect	N/A
	BGL monitor <i>Netting, top right</i>	1	11370	Consumable	N/A
	Disposable lancets <i>Stowed in BGL case</i>	5	11367	Consumable	N/A
	Test strip <i>Stowed in BG case</i>	1 x Tub	11335	Consumable	Use till expiry date
	Sharps container (0.3L) <i>Netting, bottom right</i>	1	10249	Consumable	N/A
	Thermometer <i>Netting, top middle</i>	1	11369	OET Defect	N/A
	Thermometer probe covers <i>Netting, top middle</i>	1 x box	11334	Consumable	Use till expiry date
	Viral bags <i>Netting, bottom left</i>	5	2262	Consumable	N/A

Module 1 (cont.)

	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Pen light <i>Netting, left, top</i>	1	11285	Consumable	N/A
	Fob watch <i>Netting, left, top</i>	1	11253	Consumable	N/A
	Safety specs <i>Netting, bottom middle</i>	1	3381	Consumable	N/A
	Pulse oximeter <i>Netting, top right</i>	1	11218	OET Defect	N/A
	Black pen <i>Elastic loops</i>	1	11397	Consumable	N/A
	Chinagraph pencil <i>Elastic loops</i>	1	86887	Consumable	N/A


Module 2 (Life Support)



	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Lifepak 1000 defibrillator <i>Main compartment, right side</i>	1	11354	OET Defect	N/A
	Defibrillator inner pouch <i>Main compartment, right side</i>	1	11352	OET Defect	N/A
	Defibrillator pads – adult 1 x set fitted to AED Spare set, <i>bottom main compartment, left side</i>	2	11368	Consumable	7 months before Expiry Date
	Disposable towel <i>Bottom main compartment, left side</i>	4	Station supply	Consumable	N/A

	Disposable razor <i>Lid, netting, bottom</i>	3	11366	Consumable	N/A
---	---	---	-------	------------	-----

Module 2 (cont'd)

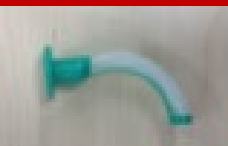	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Pocket mask – adult <i>Lid, netting, bottom</i>	1	10716	Consumable	N/A
	Pocket mask – Infant (Size 1) <i>Lid, netting, top</i>	1	11365	Consumable	Use till expiry date
	Ribbon Gauze - (iGel Tube Tie) <i>Lid, netting, Right side bottom</i>	1	12317	Consumable	N/A
	OP airway – size 1 (ISO 6.5) <i>Lid, netting, Right side top</i>	1	8684	Consumable	7 months before Expiry Date
	OP airway – size 2 (ISO 8) <i>Lid, netting, Right side top</i>	1	8685	Consumable	7 months before Expiry Date
	OP airway – size 3 (ISO 9) <i>Lid, netting, Right side top</i>	1	8686	Consumable	7 months before Expiry Date
	OP airway – size 4 (ISO 10) <i>Lid, netting, Right side top</i>	1	8691	Consumable	7 months before Expiry Date

	NP airway – 6mm (green) <i>Lid, netting, top middle</i>	1	11203	Consumable	7 months before Expiry Date
	NP airway – 7mm (orange) <i>Lid, netting, top middle</i>	1	11204	Consumable	7 months before Expiry Date
	Lubrication sachet <i>Lid, netting, top middle</i>	2	11205	Consumable	Use till expiry date

	Shears <i>Lid, netting, bottom, middle</i>	2	11252	Indent	N/A
	iGel Size 1 <i>Main compartment, left side</i>	1	12310	Consumable	7 Months before Expiry Date
	iGel Size 1.5 <i>Main compartment, left side</i>	1	12311	Consumable	7 Months before Expiry Date
	iGel Size 2 <i>Main compartment, left side</i>	1	12312	Consumable	7 Months before Expiry Date
	iGel Size 2.5 <i>Main compartment, left side</i>	1	12313	Consumable	7 Months before Expiry Date
	iGel Size 3 <i>Main compartment, left side</i>	1	12314	Consumable	7 Months before Expiry Date

	iGel Size 4 <i>Main compartment, left side</i>	1	12315	Consumable	7 Months before Expiry Date
	iGel Size 5 <i>Main compartment, left side</i>	1	12316	Consumable	7 Months before Expiry Date
	Suction <i>Main compartment, left side</i>	1	11348	Consumable	7 months before Expiry Date
	Paediatric catheter <i>Lid, netting, top</i>	1	11371	Consumable	7 months before Expiry Date
	Medical check cards (A5 sized) <i>Stored loose on top of Defib</i>	1	11403	Consumable	N/A
	Metronome <i>Clipped to defib Case</i>	1	11670	Indent	N/A
	Metronome battery CR2032	1	1030	Consumable	N/A

Module 3



	Item	Quantity	Store Number	Replacement	Replacement date
	Oxygen cylinder	1	One for one exchange	Control/SCAS (see note 6/7)	N/A
	Oxygen base pouch	1	11341	OET Defect	N/A
	Large inner pouch	2	11338	OET Defect	N/A

	Medium inner pouch	2	11339	OET Defect	N/A
	Small inner pouch	2	11340	OET Defect	N/A
	Collars (pocket behind large inner pouches)	1	6018	Consumable	N/A
	SAM splint extra large	1	11360	Consumable	N/A
	SAM splint large	1	11361	Consumable	N/A
	Triage slapper – Red P1	5	11342	Consumable	N/A
	Triage slapper – Yellow P2	5	11343	Consumable	N/A
	Triage slapper – Green P3	5	11344	Consumable	N/A
	Triage slapper – White dead	3	11345	Consumable	N/A

Large inner
pouch 1



	Item	Quantity	Store Number	Replacement	Replacement date
	Adult BVM and mask	1	10386	Consumable	7 months before Expiry Date
	Paediatric BVM and mask	1	11372	Consumable	7 months before Expiry Date
	50ML or 20ML Syringe	1	11632	Consumable	N/A
	Catheter Mount	1	12319	Consumable	7 months before Expiry Date

Large inner
pouch 2



	Item	Quantity	Store Number	Replacement	Replacement date
	Cling film + dispenser <i>*when dispenser is empty replace whole item</i>	1	10471	Consumable	N/A
	Blizzard blankets <i>*if soiled replace</i>	2	11140	Consumable	N/A
	T Pod / Pelvic Binder	1	T Pod	Consumable	N/A
	Cas Straps (2 long + 2 short) <i>*if straps are soiled post normal cleaning protocols consider replacement</i>	set	11139	Consumable	N/A

Medium inner pouch 1



	Item	Quantity	Store Number	Replacement	Replacement date
	Adhesive tape	1	11094	Consumable	N/A
	X/Large Dressing Dressing No3 (28 x 20)	2	2186	Consumable	7 months before Expiry Date
	Large Dressing *Dressing No9 (18 x 18) or Dressing No2 (20 x 15)	2	2184	Consumable	7 months before Expiry Date
	Medium dressing *Dressing No8 (12 x 12) or Dressing No1 (12 x 10)	2	2183	Consumable	7 months before Expiry Date
	Dressing No16 (eye)	2	2185	Consumable	7 months before Expiry Date

	Triangular bandage	2	2187	Consumable	7 months before Expiry Date
	Various adhesive plasters	40	2176 2177 2198 2197	Consumable	Use till expiry date
Medium inner pouch 1 (cont'd)					
	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Antiseptic wipes	10	2180	Consumable	Use till expiry date
	Tissues (small packet)	1	2181	Consumable	N/A
	Tweezers	1	2190	Consumable	N/A
	Small scissors	1	11349	Consumable	N/A

Medium inner pouch 2



	Item	Quantity	Store Number	Replacement	Replacement date
	T4 trauma dressing	1	11080	Consumable	Use till expiry date
	T6 trauma dressing	1	11081	Consumable	Use till expiry date
	OLAES modular dressing	1	11201	Consumable	No expiry date if item is still vacuum sealed
	Blast dressing	1	11202	Consumable	No expiry date if item is still vacuum sealed
	Tourniquet	2	11200	Consumable	N/A
	Chest seal	1	11346	Consumable	7 months before Expiry Date

Small inner pouch 1



	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	<i>Replacement date</i>
	Adult O2 therapy mask	1	3270	Consumable	Use till expiry date
	Child O2 therapy mask	1	9654	Consumable	Use till expiry date
	Nasal cannula	1	11364	Consumable	Use till expiry date
	Resuscitator mask	1	2195	Consumable	Use till expiry date

Small inner pouch 2



	<i>Item</i>	<i>Quantity</i>	<i>Store Number</i>	<i>Replacement</i>	
	Eye wash	1	2182	Consumable	Use till expiry date
	Glucogel	1	11353	Consumable	Use till expiry date
	Aspirin (soluble) 300mg	Pack	11401	Consumable	Use till expiry date
	Packing gauze	2	11363	Consumable	Use till expiry date

Standard Testing

Standard testing and direction on how to maintain all medical equipment is detailed within the equipment data sheets located on SharePoint.

Listed below is the general maintenance schedule for the medical response bag. Please note that the Lifepak defibrillator requires an additional **Weekly** check.

Personnel should be aware that items of medical equipment located within the medical response bag have a specific EDS.

EDS - 11/1/1 Medical Response Bag

Schedule G (Maintenance)

The bag and contents should be tested at the following intervals:

- On Acceptance
- Weekly (**Lifepak 1000**)
- Monthly
- After Use

Maintenance checks

- Carry out a full visual inspection to ensure the bag is clean and undamaged (if necessary clean with soapy water or appropriate disinfectant if contaminated).
- Check that the inventory of equipment is complete.
- Visually check each piece of equipment and if necessary, replace any missing or damaged items.
- Check that all equipment & medications (including **medical gases**) are within manufacturers use by dates. (**Consumable items with 7 months before expiry date are to be returned to stores**)
- **Blast/Oleas** dressings no longer have a use by date. The procedure is to inspect the packaging and if intact, with no loss of vacuum seal, then the dressing can roll on to the next standard test date. This process can be continued indefinitely.
- Check the O2 & Entonox cylinder as detailed in EDS/11/1/2 & EDS 11/2/4, Check Pocket Mask as detailed on EDS/11/1/4 (at 6-month intervals).

Defibrillator checks in addition to monthly maintenance.

The LP1000 performs an automatic self-test at 0300hrs daily. An alarm will sound, and the spanner symbol will be displayed if a fault is detected.

Weekly test:

- Inspect for signs of damage,
- Check 'OK' symbol is visible,
- Check battery indicator shows at least one bar,
- Check expiry dates on pads and battery,
- Check pads are pre-connected to the AED, that they are sealed and in date,
- There should be one additional set of adult pads.
- Spanner symbol in window indicates maintenance required.